

REBECCA®

Rapid Enterobacteriaceae Escherichia Coli Coliform Agar

REBECCA® is a method allowing the direct enumeration of:

- *E.coli*,
- *E.coli* & enterobacteria simultaneously
- *E.coli* & coliforms simultaneously without any confirmation



Distribuciones Industriales y Científicas S.L.

ISO16140
validation

● *E.coli*

- Direct enumeration of *E.coli* with REBECCA® base
- No confirmation needed
- ISO16140 Validation (inoculation: surface, poured plate or spiral method)
- CE2073/2005 compliant
- Results within 24h incubation at 37° C

ISO16140
validation

● *E.coli* & enterobacteria


Save money!

- Direct enumeration of *E.coli* and enterobacteria simultaneously on a single plate, by dilution of REBECCA® EB
- No confirmation needed
- ISO16140 validation (inoculation: surface, poured plate or spiral method)
- CE2073/2005 compliant

Easy to use!

- Results within 24 h incubation at single temperature of 37 °C
- Very easy enumeration of enterobacteria: **count all colonies**


* REBECCA® EB : REBECCA® base + REBECCA® EB supplement (also available in ready-to-use plates, see ordering information on the back)

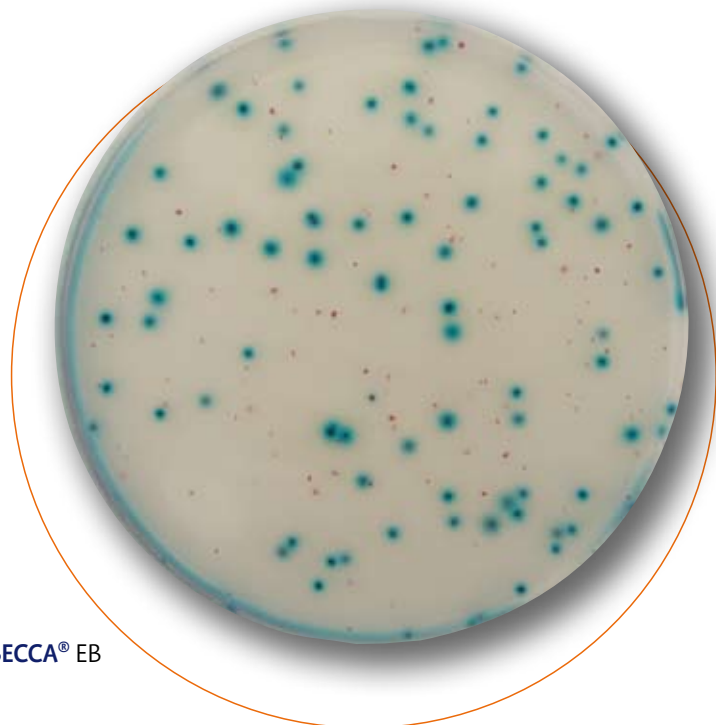

AES **CHEMUNEX**
The Rapid Microbiology Company

● Easy-to-read results

- ***E.coli***: blue colonies _____
- **Other enterobacteria**: pink to red colonies _____

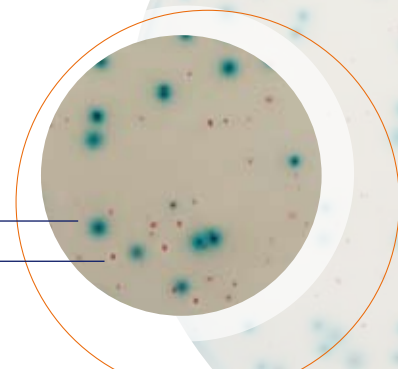


● *E.coli* & other coliforms

■ **Direct enumeration** of *E.coli* and coliforms without any confirmation

■ Results within 24h incubation at 37° C

■ **E.coli:** blue colonies

■ **Other coliform:** pink to red colonies

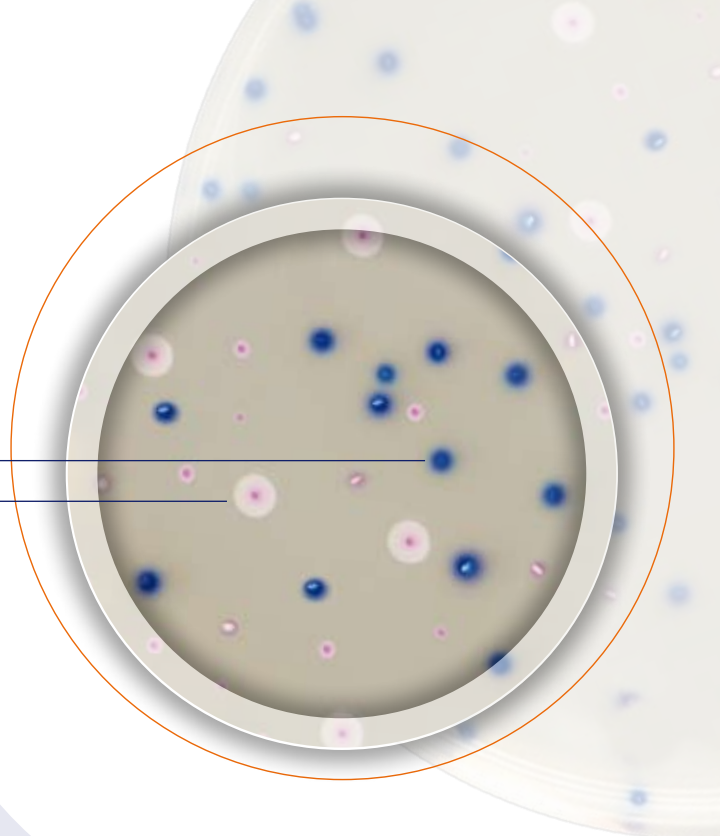

● REBECCA®

Ordering information

AEB620027 : **REBECCA®** Base
6 bottles 200 ml

AEB184135 : **REBECCA®** EB
supplement
(enterobacteria supplement)
QS for 1,2 litres

AEB520020 : **REBECCA®** EB
(ready-to-use medium for *E.coli*
and enterobacteria)
20 plates Ø 90 mm

AEB620037 : **REBECCA®** CF
(base for *E.coli* and coliforms)
6 bottles 200 ml

> Buffered peptone water

AEB610304 : 6 bottles 90 ml
AEB610308 : 6 bottles 225 ml
AEB910303 : Dilubag® 3 litres
AEB910305 : Dilubag® 5 litres
AEB140302 : 500 g

> Tryptone salt

AEB611494 : 6 bottles 90 ml
AEB611496 : 6 bottles 225 ml
AEB911493 : Dilubag® 3 litres
AEB911495 : Dilubag® 5 litres
AEB141492 : 500 g



Distribuciones Industriales y Científicas S.L.



Almeria

902.10 43 07

Fax 950 55 34 34

Ventas@dicsa.es

Cordoba

902.10 43 07

Fax 957 32 25 49

cordoba@dicsa.es

Granada

902.10 43 07

Fax 958 13 34 97

granada@dicsa.es

Jaen

902.10 43 07

Fax 953 75 83 83

jaen@dicsa.es

Malaga

902.10 43 07

Fax 952 35 63 60

malaga@dicsa.es

Sevilla

902.10 43 07

Fax 954 92 17 37

sevilla@dicsa.es